




Tracer™ DO PockeTester

Code 1761

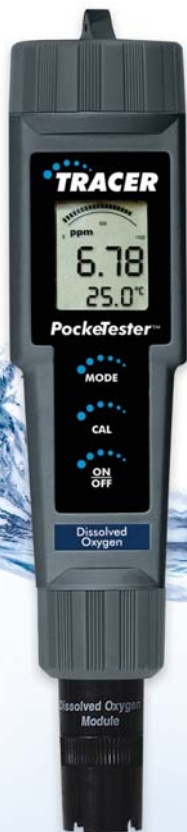

FILL - Before first use
POLARIZE - Weekly
CALIBRATE - Daily
See Pages 9-12

Water Testing Leader Since 1919!

WARNING: Reagents marked with an * are considered to be potential health hazards. To view or print a Safety Data Sheet (SDS) for these reagents go to www.lamotte.com.

Search for the four digit reagent code number listed on the reagent label, in the contents list or in the test procedures. Omit any letter that follows or precedes the four digit code number. For example, if the code is 4450WT-H, search 4450. To obtain a printed copy, contact LaMotte by email, phone or fax.

Emergency information for all LaMotte reagents is available from Chem-Tel: (US, 1-800-255-3924) (International, call collect, 813-248-0585).

WARNING! This set contains chemicals that may be harmful if misused. Read cautions on individual containers carefully. Not to be used by children except under adult supervision.

TRACER™
DISSOLVED OXYGEN • CODE 1761

Introduction	4
Specifications	4
Contents	4
Parts & Accessories	5
Meter Description	
Front Panel Description	5
Display	5
Operation	
Filling the Electrode	6
Powering the TRACER	6
Startup Polarization Period	6
Turn-on Diagnostics	6
Changing the Temperature Units	6
Changing Salinity Compensation	7
Changing the Altitude Compensation	7
Changing Measurement Units	7
Auto-Power Off Feature	8
Disabling Auto-Power Off Feature	8
Low Battery Indicator	8
Storing Readings	8
Recalling Stored Readings	8
Clearing Stored Readings	8
Set Up	
Membrane Cap Assembly - Filling and Replacement	9
Polarization	10
Calibration	11
Optional Zero Oxygen Calibration	11
Dissolved Oxygen Measurement	12
Operational Matrix	13
Troubleshooting	16
Maintenance	
Battery Replacement	17
Battery Safety	17
Electrode Replacement	17
Replacing the Membrane Cap	17
Electrode Storage	18
Warranty	18
Short Form Instructions (back cover)	

INTRODUCTION

Congratulations on your purchase of the Dissolved Oxygen/Temperature TRACER. Units of measure are % saturation, mg/L or ppm for dissolved oxygen and C° or F° for temperature. Advanced features include data hold, 25 point memory, auto shut off, automatic temperature compensation, and salinity/altitude compensation. Careful use and maintenance will provide years of reliable service.

SPECIFICATIONS

<i>Display</i>	<i>2000 count, Dual function 3 ½ digit LCD with Bar graph</i>
<i>Display size</i>	<i>24 mm x 20 mm</i>
<i>Sensor</i>	<i>Polarographic type</i>
<i>Membrane</i>	<i>Bonded membrane cap with threaded fitting</i>
<i>Operating Temp.</i>	<i>Range 0 to 50°C (32 to 122°F)</i>
<i>ATC Range</i>	<i>0 to 50°C (32 to 122°F)</i>
<i>Salinity Compensation</i>	<i>0 to 50 ppt in 1 ppt (part per thousand) increments</i>
<i>Altitude compensation</i>	<i>6000 m (0 to 20,000 ft) in 300m (1000 ft) increments</i>
<i>Measurement Storage</i>	<i>25 tagged (numbered) data sets with recall</i>
<i>Battery Power</i>	<i>Four (4) CR2032 button batteries</i>
<i>Low Battery Indication</i>	<i>'BAT' appears on the LCD</i>
<i>Auto Power Off</i>	<i>After 10 minutes of inactivity (APO override available)</i>
<i>Dimensions/Weight</i>	<i>36 x 173 x 41mm (1.4 x 6.8 x 1.6"); 110g (3.8 oz)</i>

Measurement	Range	Resolution	Accuracy
<i>% Saturation</i>	<i>0 to 200.0%</i>	<i>0.1%</i>	<i>±2.0% FS (full scale)</i>
<i>Dissolved Oxygen Concentration</i>	<i>0 to 20.00 mg/l</i>	<i>0.01 mg/l</i>	<i>±2% FS</i>
<i>Temperature</i>	<i>0 to 50 °C</i>	<i>0.1 °C</i>	<i>±1.0 °C</i>
	<i>32 to 122°F</i>	<i>0.1 °F (0 to 99°F);</i>	<i>±1.8°F</i>
		<i>1.0 °F (>100°F)</i>	

CONTENTS

Dissolved Oxygen TRACER Kit *Code 1761*

Includes:

- DO TRACER Body*
- DO Bonded Membrane Cap Assembly*
- Electrolyte Solution*
- Sample cup and Lanyard*

PARTS & ACCESSORIES

<i>DO Replacement Sensor Kit (sensor module, membrane cap assembly, electrolyte solution)</i>	<i>Code 1762</i>
<i>DO Membrane Cap Kit (6 membrane cap assemblies, electrolyte solution, Polishing Strips)</i>	<i>Code 1761M</i>
<i>DO Extension Cable, 1 meter (with probe guard and weight)</i>	<i>Code 1763</i>
<i>DO Extension Cable, 5 meter (with probe guard and weight)</i>	<i>Code 1764</i>


METER DESCRIPTION

Front Panel Description

1. Battery compartment cover
2. LCD Display
3. MODE button
4. CAL button
5. ON/OFF button
6. Electrode Retaining Collar
7. Dissolved Oxygen Sensor
8. Bonded membrane cap assembly
9. Membrane & cathode
(Note: The Electrode storage cap is not shown)




Bottom View of Electrode

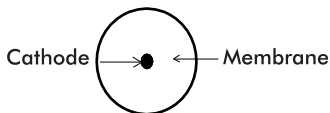

Display

1. Bar graph display
2. Measurement units
3. Main display
4. Low Battery indicator
5. Temperature display
6. Reading Hold indicator
7. Polarization indicator (meter off)

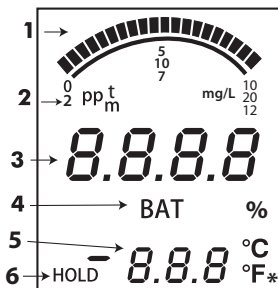

OPERATION

Filling the Electrode

The electrode is shipped dry and requires filling with the supplied electrolyte solution before the first use. The membrane should be in place and does not need to be replaced. Follow the procedure on page 9 for filling the membrane cap.

Powering the Tracer

The Tracer uses four (4) CR2032 Lithium Ion Batteries. Press the **ON/OFF** button to turn the meter on or off. If the batteries are weak, the 'BAT' indicator will appear on the display. The auto power off feature shuts the Tracer off automatically after approximately 10 minutes of inactivity. The auto power off feature may be temporarily disabled for convenience or for extended polarization time.

Startup Polarization Period

When the DO Tracer is powered up for the first time, the electrode requires polarization. In order for this to occur a polarization period of approximately three minutes must elapse before measurements can be made. After the meter has been turned on, a special circuit in the meter will maintain a very small biasing current to the electrode for a period of seven days. This will keep the electrode polarized and enables the user to make immediate measurements within a 7 day period without waiting for the electrode to re-polarize. Each time the Tracer is turned on, the polarization timer is reset and the 7 day polarization period is initiated. A small asterisk at the bottom right corner of the display will indicate that the polarization timer circuit is active only when the meter is off. If the Tracer is not used for more than seven days, a 3 minute waiting period is required before a measurement can be made.

Turn-On Diagnostics

1. When the meter is turned on 'SELF' and 'CAL' will be displayed while the meter runs a diagnostic routine.
2. During this time the meter is recalling the User Calibration data, performing self diagnostics & initializing the circuitry.
3. After this function is complete, the meter proceeds to the normal measurement mode.

Changing the Temperature Units

To change the displayed temperature units between °F and °C

1. With the unit OFF, press and hold **CAL**.
2. With **CAL** depressed, momentarily press **ON/OFF** to turn the Tracer ON.
3. Release **CAL**. 'SELF CAL' will be displayed.
4. Repeat Steps 1 – 3 to change to the other temperature unit.

Changing the Salinity Compensation

1. With the Tracer ON, momentarily press **CAL** twice in quick succession. 'SAL' will be shown on the lower part of the display.
2. Momentary press **MODE**. Each press of **MODE** increases the salinity compensation by 1ppt (part per thousand); the available range is 0 to 50ppt.
3. Momentarily press **CAL** to save the compensation setting and return to the normal measurement mode.

Changing the Altitude Compensation

1. With the Tracer ON, momentarily press **CAL** twice in quick succession. 'SAL' will be shown in the lower part of the display.
2. Press and hold **CAL** again for 2 seconds to enter the Altitude Mode. 'Ald' will be shown in the lower part of the display.
3. The factory default value is sea level. Each press of **MODE** increases the altitude compensation by 1000ft. The maximum value is 20 presses (20,000ft above sea level).
4. Momentarily press **CAL** to save the compensation setting and return to the normal measurement mode.

Changing the Measurement Units

The meter can be set to measure % saturation, dissolved oxygen in mg/l, and dissolved oxygen in parts per million (ppm).

To change the mode:

1. Press and hold **MODE** for 2 seconds and the display will begin to scroll through the available units of measure :
 - % saturation
 - dissolved oxygen, mg/L (milligrams per liter)
 - dissolved oxygen, ppm (parts per million)
2. When the desired units are displayed, release **MODE** and the unit will return to normal operating mode.

Note: The 'HOLD' function cannot be on when changing the measurement function. If 'HOLD' is displayed in the lower left corner of the display, briefly press the **MODE** button to turn it off.

Auto-Power Off Feature

The auto power off feature automatically shuts the meter off 10 minutes after the most recent button press. To disable this feature, refer to the Disable Auto-Power Off section.

Disable Auto-Power Off Feature

With the meter **ON**, momentarily press **CAL**, then quickly press and hold both **MODE** and **ON/OFF** until 'oFF' is displayed. To restore the Auto Power Off Feature turn the meter off and on again using **ON/OFF**.

Low Battery Indicator

When the battery voltage falls below the operating threshold, 'BAT' will appear on the display. Refer to the Maintenance section for battery replacement information.

Storing Readings

1. Press **MODE** to store a reading. The storage location number will be shown on the lower display and the main display will show the stored reading. The meter will enter the HOLD mode and 'HOLD' will appear on the display.
2. Press **MODE** again to exit the HOLD mode and return to normal operation. The next time **MODE** is momentarily pressed another reading will be stored and so on.
3. If more than 25 readings are stored, previously stored readings, starting with reading number 1, will be overwritten.

Recalling Stored Readings

1. Momentarily press the **CAL** and then, within 4 seconds, momentarily press **MODE**. The last stored data point location will be displayed (1 to 25). Each time **MODE** is momentarily pressed the next most recently stored data point will be displayed.
2. After the last data point stored is displayed, press **MODE** again to return the display to the beginning of the list.
3. Press **CAL** at any time to stop the data retrieval process and return the meter to the normal measurement mode.

Clearing Stored Readings

1. With the meter **ON** press and hold **ON/OFF** for 4 seconds.
2. When 'clr' is shown on the display the memory is cleared.

SET UP


Membrane Cap Assembly - Filling and Replacement

The electrode is shipped dry and requires filling with the supplied electrolyte solution before the first use. The membrane should be in place and does not need to be replaced.

Important Note: Do not touch the membrane as skin oils will interfere with the oxygen permeability rate of the membrane.

1. The electrode should remain attached to the meter during the process.
2. Carefully remove the membrane cap assembly from the electrode by firmly unscrewing it in a counter-clockwise direction.
3. Discard the used membrane cap assembly.

NOTE: When a membrane cap is installed the membrane must tightly stretch over the cathode. If a membrane cap is removed, it can not be reinstalled because the membrane will no longer stretch properly over the cathode.

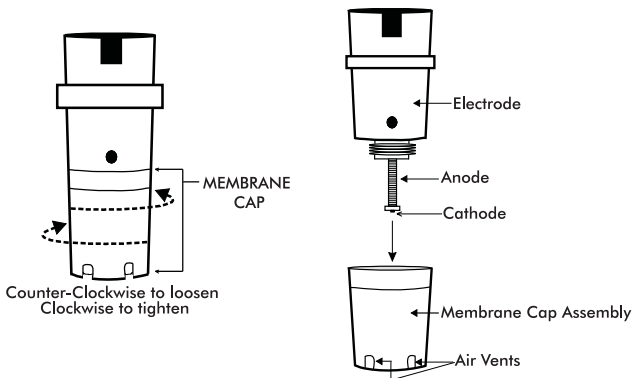

4. Rinse the cathode and anode with deionized water.
5. Set the membrane cap assembly on a flat surface.
6. Fill the membrane cap assembly with Electrolyte Solution (DO600-EL) to the bottom of the threads on the inside of the cap. To remove air bubbles hold the cap steady with one hand and tap the side of the assembly sharply with a pen.

7. Do not move the membrane cap assembly. Dip the electrode into the assembly and remove it. Repeat this step several times. With each dip, progressively push the electrode deeper into the electrode solution in the assembly. This dipping technique will minimize the introduction of air bubbles into the electrolyte solution. Air bubbles can compromise dissolved oxygen measurements.
8. Carefully pick up the membrane cap assembly and slowly screw it onto the electrode in a counter-clockwise direction until it is fully tightened. While tightening the cap, excess electrolyte solution will leak out. This is normal and desirable since it will minimize air pockets in the solution. Air pockets can compromise dissolved oxygen measurements. If the electrode solution does not overflow, an insufficient amount was added. Repeat the process with a new membrane cap assembly if air bubbles or air pockets are visible when the meter is inverted. Membrane cap assemblies can not be reused.

NOTE: Air bubbles outside the cathode (6 mm in diameter in the center of the probe) will not interfere.

9. Rinse excess electrolyte from the outside of the membrane cap assembly and the electrode.
The pH of the electrode filling solution in polarographic dissolved oxygen meters changes with time resulting in less accurate results. It is recommended that the membrane and electrode filling solution be changed every 2 - 4 weeks depending on usage. Because of the configuration of the membrane and electrode, the membrane is stretched by the cathode and cannot be reused. It must be replaced.

Polarization

If the meter has not been used for seven (7) days or longer the electrode will not be fully polarized. A fully polarized electrode is indicated by a small, flashing, asterisk in the lower right corner of the display when the meter is OFF. If the asterisk is absent, turn the meter ON and allow the electrode to fully polarize. This may take 2-3 minutes. The asterisk will appear to indicate full polarization.

Calibration

Calibration should be performed daily.

(Be sure the electrode is filled and fully polarized before calibrating)

1. Turn the meter ON.
2. Press and hold **MODE** until '%' is displayed.
3. Remove the electrode cap. The electrode membrane must be clean and dry or the calibration will be incorrect. Do not touch the membrane. Skin oil will affect the electrode response.
4. Moisten the sponge in the electrode cap with distilled water or tap water. Do not soak the sponge.
5. Replace the electrode cap on the electrode. There must be an air gap between the membrane and the sponge.
6. Wait until the reading stabilizes. Press and hold **CAL** until 'CAL' is shown on the display. The readings will blink '101.7' and 'SA' will appear. 'SA' will not appear if the calibration fails.
7. When the calibration is complete 'End' will appear and the meter will return to the measurement mode.

Optional Zero Oxygen Calibration

The zero oxygen calibration will improve the measurement accuracy for very low or very high dissolved oxygen samples.

1. Place the electrode in a zero dissolved oxygen calibration solution, such as 5% sodium sulfite.
2. Slowly stir the sample with the electrode and wait for the reading to stabilize. This may take a while, depending on the electrode history.
3. Press **CAL** until 'CAL' is shown on the display. The readings will blink '101.7' and 'SA' will appear. 'SA' will not appear if the calibration fails.
4. When the calibration is complete 'End' will appear and the meter will return to the measurement mode.

NOTE: Sodium Sulfite can become deposited on the electrode and on the coined surface of the electrode retaining collar. The presence of the Sodium Sulfite will negatively affect dissolved oxygen measurements unless it is removed completely from the electrode.

DISSOLVED OXYGEN MEASUREMENT

(Be sure the electrode is filled and fully polarized. Calibrate daily.)

1. Cover the electrode with the electrode cap. The sponge in the electrode cap should be moistened (not soaked) with distilled water or tap water.
2. Press **ON/OFF** to turn the meter ON. The meter will self calibrate.
3. Press and hold **MODE** to select the desired unit of measure.
4. Remove the electrode cap.
5. Place the electrode in the sample to be measured. Stir the sample with the electrode.
6. Wait for the reading to stabilize.
7. Press **ON/OFF** to turn the meter OFF. Rinse the electrode with distilled water or tap water. Replace the electrode cap.

NOTE:

- The unit of measure selected when the meter was last turned OFF will be retained.
- Stir the sample with the electrode to remove any air bubbles that may be trapped on the membrane surface and to obtain accurate results.
- Do not submerge the electrode above the bottom edge of the electrode collar.
- For maximum accuracy, allow sufficient time for the temperature of the probe to reach the temperature of the sample before taking a reading. The larger the difference in temperature between the electrode and the sample the longer it will take for the reading to stabilize. Stabilization will be indicated by a stable temperature reading on the display. Stabilization time can vary from ten (10) seconds to five (5) minutes.
- The DO Tracer uses a polarographic sensor, which consumes oxygen at the sensing surface. This requires a constant movement of the sample across the membrane to maintain a constant dissolved oxygen level. It is recommended that the sample be stirred gently with the electrode as the measurement is made, or if in a laboratory environment, the sample should be stirred using a stir plate.

OPERATIONAL MATRIX

Function / Resulting Action	Power Status	Mode Setting	Required Button Press Sequence	Comments
ON/OFF	Any	Any	Momentary press of the ON/OFF button	
Water Saturated Air Calibration	On	Any	Place electrode in electrode cap. Press & hold CAL button for 2 seconds.	
Zero Calibration	On	Any	Place electrode in zero DO solution. Wait for stability. Press & hold CAL button for 2 seconds.	Also works with probe removed. (Cal at zero current)
Store Reading	On	Any	Momentary press of MODE button.	Stores and holds reading. 'HOLD' displayed
Hold Release	On	While In Hold Mode	Momentary press of the MODE button.	
Enter Memory Retrieval	On	Any	Momentary press of CAL button followed by a momentary press of the MODE button within 4 seconds.	If no data is stored in memory, 'End' is displayed briefly, and then returns to last mode.
Scroll Stored Readings	On	Memory Recall	Momentary press of MODE button.	
Exit Memory Retrieval	On	Memory Recall	Momentary press of CAL button.	
Clear Stored Memory	On	Any Measurement Mode	Press and hold ON/OFF button for 4 seconds	'clr' is displayed.

Function / Resulting Action	Power Status	Mode Setting	Required Button Press Sequence	Comments
Change Measurement Mode	On	Any	Press and hold MODE button for at least two seconds. Modes will scroll until button is released.	
Enter Salinity Compensation	On	Any	Press and release the CAL button twice quickly in succession.	'SAL' is displayed.
Changing Salinity Compensation	On	SAL	Momentary press of MODE button.	Each button press increases the ratio by 1 ppt (part per thousand). The value cycles from 0 to 50 ppt.
Exit Salinity Compensation	On	SAL	Press and release CAL button for 2 seconds to enter Altitude Compensation. Or press and release CAL button one more time to enter measurement mode.	It is necessary to press CAL button to save changes. If the unit times out, no changes are saved.
Enter Altitude Compensation	On	Any or SAL	Press CAL button twice quickly. Unit enters Salinity Mode. Press CAL button for 2 seconds to enter Altitude Compensation Mode.	'Ald" is displayed. Unit times out in 5 seconds if no button is pressed, reverts to previous mode.

Function / Resulting Action	Power Status	Mode Setting	Required Button Press Sequence	Comments
Changing Altitude Compensation	On	Ald	Momentary press of the MODE button.	Each button-press increases altitude by 1,000 ft. The value cycles from 0 to 20,000 ft.
Exit Altitude Compensation	On	Ald	Momentary press CAL button to exit and save changes.	It is necessary to press the CAL button to save changes. If the unit times out, no changes are saved.
Change Temperature Units	Off	OFF mode	Press and hold the CAL button then momentarily press the ON/OFF button. Release the CAL button after 'SELF CAL' appears.	
Override Auto Power Off	On	Any	Press CAL button momentarily then simultaneously press and hold MODE and ON/OFF buttons for 2 seconds.	
Default Reset	Off	OFF mode	Simultaneously press ON/OFF , CAL and MODE buttons momentarily.	'dFLt' will be displayed.

TROUBLESHOOTING


- If the display appears to be frozen, it is possible that the Data Hold mode has been inadvertently activated. This is indicated by 'HOLD' in the bottom left of the display. To exit the Data Hold mode press **MODE** or turn the meter off and back on again.
- If the meter locks up and pressing any buttons fails to revive it, remove the batteries and then reinsert the batteries.

Symptom	Possible Cause	Action
Unit will not power on	Batteries not in place Dead batteries Batteries polarity incorrect	Insert, replace, or reorient batteries
“BAT” indicator shown on display	Batteries are weak	Replace batteries
Unstable readings	Insufficient electrolyte in probe. Air bubbles visible when probe inverted. Electrolyte is depleted	Replace electrolyte, and membrane cap assembly
Readings drift down	Insufficient stirring	Stir sample with meter or move probe in sample
Slow response	Dirty or damaged membrane	Replace electrolyte and membrane cap
Electrode can not be calibrated	Depleted electrolyte Dirty or damaged membrane	Replace electrolyte and membrane cap
Electrode can not be calibrated after replacing electrolyte and membrane cap	Dirty probe. Cathode is not a shiny gold color.	Clean cathode with cleaning cloth or cleaning paper. Do not over polish.
Sample reading is frozen	Meter is in HOLD mode Meter is locked	Release HOLD by momentarily pressing MODE button. Remove batteries. Press ON/OFF button. Replace batteries and press ON/OFF button. Stored readings will be lost.

MAINTENANCE

Battery Replacement

1. Twist off the battery compartment cover.
2. Holding the battery housing in place with a finger, pull out the battery carrier using the two small tabs.
3. Replace the four (4) CR2032 batteries observing proper polarity.
4. Replace the battery carrier, reattach the battery compartment cap and tighten securely.

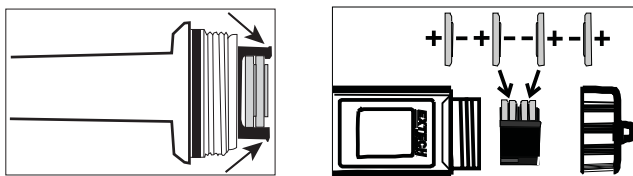

Battery Safety

- Dispose of batteries responsibly; observe local, state, and federal regulations.
- Never dispose of batteries in a fire; batteries may explode or leak.
- Never mix battery types; install new batteries of the same type.

Electrode Replacement

1. Turn the meter OFF. Turn the electrode retaining collar clock-wise and remove it.
2. Gently rock the electrode from side to side, pulling it away from the meter until it disconnects.
3. Align the positioning “keys” of the replacement electrode and the main body housing. Carefully push the electrode into the meter socket until it is fully seated. CAUTION: Take care to align pins carefully. Bent or broken pins will cause the meter to malfunction.
4. Tighten the electrode retaining collar until the electrode is firmly sealed to the meter.

Replacing the Membrane Cap Assembly

The pH of the electrode filling solution in polarographic dissolved oxygen meters changes with time resulting in less accurate results. It is recommended that the membrane and electrode filling solution be changed every 2 - 4 weeks depending on usage. Because of the configuration of the membrane and electrode, the membrane is stretched by the cathode and cannot be reused. It must be replaced. See page 9.

Electrode Storage

Store the TRACER with the Electrode Storage Cap in place. The sponge in the cap should always be moist without excess liquid.

Electrode Cleaning Recommendations

Contaminant	Cleaning Solution	Procedure
Salt build-up	Deionized water	Soak in 20% or weaker Acetic Acid solution for 1 minute. Rinse thoroughly with DI water.
Oil	Warm water and household detergent	Clean with mild detergent solution. Rinse thoroughly with DI water.

NOTE: Do not contact the membrane with the cleaning solution. Install a new membrane after cleaning the electrode.

WARRANTY

LaMotte Company warrants this instrument to be free of defects in parts and workmanship for 1 year from the date of shipment and warrants this probe to be free of defects in parts and workmanship for 6 months from the date of shipment. If it should become necessary to return the instrument for service during or beyond the warranty period, contact our Technical Service Department at 1-800-344-3100 or tech@lamotte.com for a return authorization number or visit www.lamotte.com for troubleshooting help. The sender is responsible for shipping charges, freight, insurance and proper packaging to prevent damage in transit. This warranty does not apply to defects resulting from action of the user such as misuse, improper wiring, operation outside of specification, improper maintenance or repair, or unauthorized modification. LaMotte Company specifically disclaims any implied warranties or merchantability or fitness for a specific purpose and will not be liable for any direct, indirect, incidental or consequential damages. LaMotte Company's total liability is limited to repair or replacement of the product. The warranty set forth above is inclusive and no other warranty, whether written or oral, is expressed or implied.

To register your meter with the LaMotte Service Department, go to www.lamotte.com and choose SUPPORT on the top navigation bar.

Short Form Instructions

Fill the membrane cap assembly before the first use.

Follow Set Up instructions before the first use.

Read the entire manual before the first use.

Calibrate daily.

Polarize weekly.

Dissolved Oxygen Measurement Procedure

1. Press **ON/OFF** button to turn meter ON.
2. Place TRACER in sample. Stir.
3. Wait for reading to stabilize.
4. Record result.
5. Press **ON/OFF** to turn meter off.
6. Rinse electrode.
7. Replace electrode cap.

LaMOTTE COMPANY

Helping People Solve Analytical Challenges

PO Box 329 • Chestertown • Maryland • 21620 • USA
800-344-3100 • 410-778-3100 (Outside USA) • Fax 410-778-6394

Visit us on the web at www.lamotte.com