



WIDE RANGE pH TEST KIT

OCTA-SLIDE 2, pH 5.0-10.0

CODE 3353-01

QUANTITY	CONTENTS	CODE
30 mL	*Wide Range pH Indicator	*2218-G
2	Test Tubes, 2.5-5-10 mL, plastic, w/caps	0106
1	Wide Range pH Octa-Slide 2 Bar, 5.0-10.0	3483-01
1	Octa-Slide 2 Viewer	1101


*WARNING: Reagents marked with an * are considered to be potential health hazards. To view or print a Safety Data Sheet (SDS) for these reagents go to www.lamotte.com. Search for the four digit reagent code number listed on the reagent label, in the contents list or in the test procedures. Omit any letter that follows or precedes the four digit code number. For example, if the code is 4450WT-H, search 4450. To obtain a printed copy, contact LaMotte by email, phone or fax.

Emergency information for all LaMotte reagents is available from Chem-Tel:
(US, 1-800-255-3924)
(International, call collect, 813-248-0585).

Warning! This set contains chemicals that may be harmful if misused. Read cautions on individual containers carefully. Not to be used by children except under adult supervision.

To order individual reagents or test kit components, use the specified code number.


USE OF THE OCTA-SLIDE 2 VIEWER

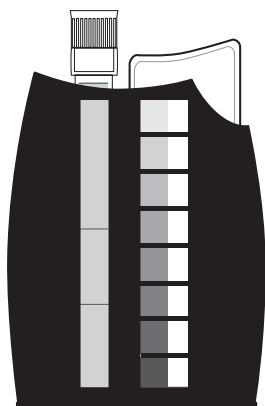

The Octa-Slide 2 Viewer should be held so non-direct light enters through the back of the Viewer. Insert the Octa-Slide 2 Bar into the Viewer. Insert the reacted sample into the top of the Viewer. Match the color of the reaction to the color standards.


PROCEDURE

1

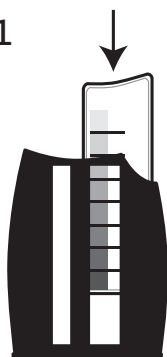


Insert Wide Range
pH Octa-Slide 2
Bar (3483-01) into
the Octa-Slide 2
Viewer (1101)

2

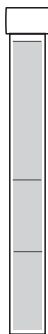


Fill a test tube
(0106) to the
10 mL line with
sample water.

3

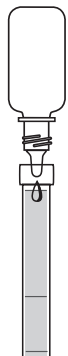


Add 8 drops of
*Wide Range pH
Indicator (2218).

4

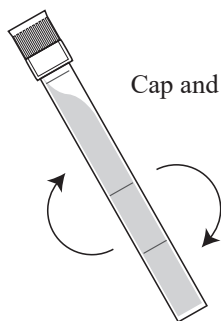


Cap and mix.

5

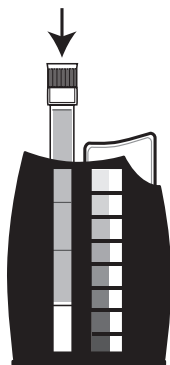

Insert test tube
into Octa-Slide 2
Viewer (1101).

6

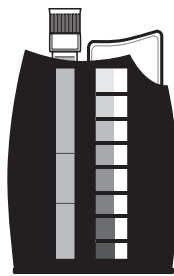

Match sample color
to a color standard.
Record as pH.

LaMOTTE COMPANY

Helping People Solve Analytical Challenges

PO Box 329 · Chestertown · Maryland · 21620 · USA
800-344-3100 · 410-778-3100 [Outside U.S.A.] · Fax: 410-778-6394
Visit us on the web at www.lamotte.com